






Rovos R653L Foot Massager



Product Operation Guide

- 
ON/OFF
 Long press to turn ON.
 Long press again to turn OFF
 Default timer setting is 15 minutes.
- 
HEAT
 Press to turn ON.
 Press again to turn OFF. 3 Heating intensity
- 
INTENSITY
 Press to switch massage INTENSITY
 3 intensity
- 
AUTO AUTO
 Press to enable and switch AUTO Mode
 1. A combination soft massage
 2. A combination moderate massage
 3. A combination strong massage

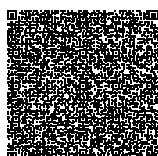


Product Parts




This product uses only the 220-240V, 50-60Hz power supply.
Make sure your indoor voltage fits this specification parameter.

- Connect the AC power cord to the socket.

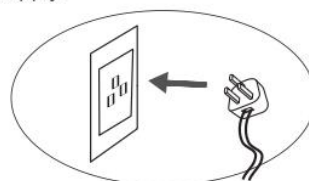

Do not insert or remove the plug with a wet hand

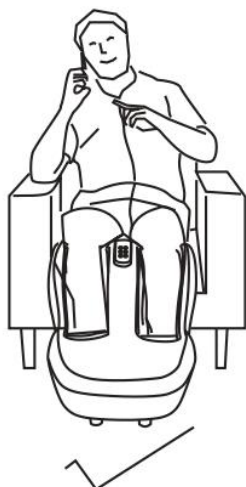

Do not put your hands inside the cloth cover.



Product Specifications

| | |
|-------------------|-----------------------|
| Power supply | AC 220V-240V, 50-60Hz |
| Power consumption | 45W |
| Timer | 15 minutes |
| Dimension (LxWxH) | 420*455(mm) |



Manufactured By:
Fujian Rovos Fitness Co., Ltd
XibeiYang Industrial Zone, Fuan, Fujian, China

Registrasi Kemenkes RI
AKL

Distributed By:
PT. Rovos Fit Indonesia
Jl. Nusantara Timur Blok D41 Sunter Sacna
Jakarta Utara Indonesia 14350

